

Whitley M45 Workstation

OVERVIEW

SPECIFICATIONS

OPTIONS

CONSUMABLES

IMAGES/VIDEOS

DOWNLOADS



Variable atmosphere and an instant access port

- Optional removable front (with instant access or traditional sleeved ports) to enable the transfer of larger batches of samples or pieces of equipment for use inside the chamber.
- Built-in rapid airlock allows the transfer of plates in just 60 seconds.
- Colour, touch-screen control panel for ease of use and for visual display of parameters such as temperature, humidity, and airlock cycle status.
- Reminder to perform calibration can be set for between 7 and 21 days.
- Save time with the automatic commissioning cycle.

The M45 variable atmosphere workstation is ideal for the study and isolation of *Campylobacter* spp, *Helicobacter pylori* and other similarly fastidious microaerophiles. The M45 has three ports for convenient user access to the entire incubation and working areas.

This is a 4-gas system with built in oxygen, CO₂ and hydrogen sensors allowing users to programme precise gas concentrations, perfect for manipulating samples in a sustainable microaerobic environment. Microaerobic workstations provide the most flexibility of any models in the Don Whitley Scientific range; allowing Anoxic/Anaerobic, Microaerobic or Hypoxic atmospheres to be created if desired without the need to change incoming gasses.

A range of optional accessories, such as automatic humidification, are available to tailor this workstation to your particular requirements.

This workstation can also be specified with the Whitley HEPA Filtration System and can accommodate up to 900 x 90mm Petri dishes..

[VIEW SPECIFICATION](#)





If a larger airlock is needed, try the M85 with its 30 litre airlock ...

WHITLEY M85 WORKSTATION

000 «Диаэм»

Москва
ул. Магаданская, д. 7, к. 3 ■ тел./факс: (495) 745-0508 ■ sales@dia-m.ru

www.dia-m.ru

С.-Петербург
+7 (812) 372-6040
spb@dia-m.ru

Новосибирск
+7(383) 328-0048
nsk@dia-m.ru

Воронеж
+7 (473) 232-4412
vrn@dia-m.ru

Йошкар-Ола
+7 (927) 880-3676
nba@dia-m.ru

Красноярск
+7(923) 303-0152
krsk@dia-m.ru


Казань
+7(843) 210-2080
kazan@dia-m.ru

Ростов-на-Дону
+7 (863) 303-5500
rnd@dia-m.ru

Екатеринбург
+7 (912) 658-7606
ekb@dia-m.ru

Кемерово
+7 (923) 158-6753
kemerovo@dia-m.ru

Армения
+7 (094) 01-0173
armenia@dia-m.ru

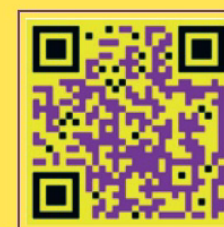

[Terms and Conditions](#)

[Cookie Policy](#)

[Privacy Policy](#)

[Sitemap](#)

sales@dwscientific.co.uk

[+44\(0\)1274 595728](tel:+44(0)1274595728)

Join our Mailing List:

JOIN MAILING LIST

Don Whitley Scientific Limited. Registered in England no. 1342672. Registered Office: Victoria Works, Victoria Street, Bingley, West Yorkshire, BD16 2NH, UK

© 2021 Don Whitley Scientific Limited Copyright: All rights reserved - No content can be reproduced without our prior written consent.

[Web Design - Rejuvenate Digital](#)